



KAROLINSKA INSTITUTET
- a medical university -

EXHIBIT 1

Department of: CLINICAL IMMUNOLOGY

Laboratory Journal No: .

Name: ~~XXXXXXXXXX~~
HANS GRENBLUND

Group:

Date: [REDACTED]
from to



cont. from page no:

Study

cont. on page no:

Project no.

Study no.

27

Assembly of a synthetic gene coding for
Tel d1 Chain 1 using Tag

PCR

Oligos 132, 133, 134, 135 10 μ M

1 μ l of each oligo 132-135

1 μ l dNTP 10 mM

0.5 μ l Tag

1 μ l 10x Tag-buffer

3.5 μ l H₂O

10 μ l

→ PCR Eppendorf program HAN51

94°C 1 min


annealing 65-65°C 1.5 min (grad)

elongation 68°C 20 min

30 cycles ending 10 min elongation

+ 4°C

Result A strong band around 300 (exp. ~ 200)

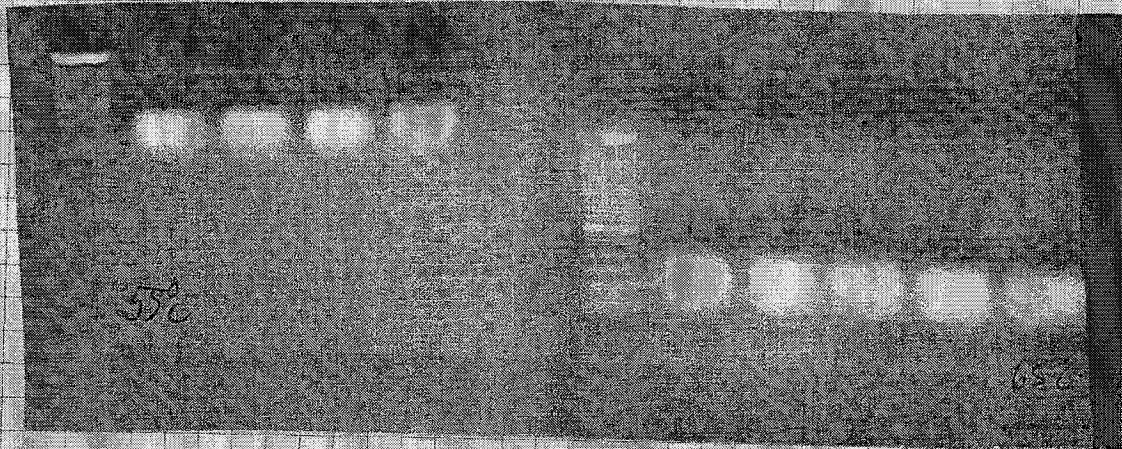

Read and understood by

Signature

Date

Signature

Date



Laboratory Journal

248910

cont. from page no.

Study

Assembly of chain 2 Fel d1

cont. on page no.

Project no.

Study no.

primers 127-131 + 138

Using 3 different DNA polymerases

Method as: assembly chain 1"

Tag
pfu

Ampli Tag

expected band
at 294



to 248913
Lane 1

Read and understood by

Signature

Date

Signature

Date



cont. from page no:

Study

Expression and purification of
Feld1 chain 1 and chain 2.

cont. on page no:

Project no.

Study no.

Feld1 chain 1 (clone 42) and Feld1 chain 2 clone 29 was ligated into pET 20b and electroporated into BL-21 DE3 pLysS after having been cut from pT7-Blue containing the correct sequence. (see binder HG1, Feld1) Sequencing of pET 42/pET 29 was done according to standard protocol (ABI) and the results can be seen on the opposite side.

Both Feld1 chain 1 (Fd1:1) and Feld1 chain 2 (Fd1:2) was expressed according to standard protocol and purified on a HiTrap (chelex) column loaded with NiSO₄.

Chain 1 was soluble, after ultra sonication chain 1 was found in 20 mM Tris-HCl pH fraction, while chain 2 was found in the inclusion bodies after "washing" with 2M urea buffer + 20 mM Tris-HCl pH 8.0. The inclusion bodies were solubilized in 6M guanidine, transferred to 6M urea buffer (20 mM Tris-HCl pH 8.0 + 0.5M NaCl) via 6M Tris-HCl trap. Purification was done on FPLC

Read and understood by

Signature

Date

Signature

Date

test serie/analyse Biedt Algeet II

| Personnummer | knr | g.1 | et prick | |
|--------------|-----|-------|----------|--------|
| 340217-3235 | 29 | 0.47 | | Amphar |
| 370413-3218 | 59 | 11.4 | 3 | Amphar |
| 380620-1236 | 50 | 8.2 | 3 | |
| 419317-3210 | 95 | 1.24 | | |
| 429508-3215 | 102 | 2.21 | | |
| 471013-3254 | 110 | 0.43 | | |
| 430308-3277 | 118 | 1.3 | | |
| 441118-3272 | 135 | 1.03 | 2 | |
| 471031-3273 | 172 | 1.07 | | |
| 480710-3218 | 193 | 0.82 | | |
| 500228-3217 | 203 | | 2 | |
| 500615-3233 | 234 | 0.73 | | |
| 510414-3231 | 235 | | 2 | |
| 511020-3212 | 250 | 0.84 | 3 | |
| 531018-3232 | 243 | 0.37 | | |
| 531108-3235 | 245 | 0.85 | 2 | |
| 540209-3214 | 253 | | 2 | |
| 540823-3232 | 251 | 87.7 | | |
| 550521-3213 | 272 | 26.5 | 2 | |
| 560701-3218 | 283 | 0.45 | | Amphar |
| 570415-3211 | 285 | 0.45 | | |
| 580611-3215 | 317 | 0.82 | 2 | |
| 600623-3215 | 341 | | 2 | |
| 640316-3216 | 383 | 2.06 | | |
| 621110-3255 | 378 | | 2 | |
| 650127-3212 | 380 | | 2 | |
| 630703-3230 | 388 | 1.56 | | |
| 640724-3234 | 398 | 0.47 | | |
| 660628-3210 | 403 | | 2 | |
| 680828-3212 | 406 | 0.76 | | |
| 680703-3241 | 416 | 3.06 | 2 | |
| 671004-3256 | 422 | 26.5 | 3 | Amphar |
| 680220-3213 | 426 | | 2 | |
| 680520-3210 | 428 | 0.50 | | |
| 690108-3235 | 434 | 21.45 | 3 | Amphar |
| 730812-3251 | 452 | 5.33 | | |
| 740300-0044 | 454 | 8.33 | 2 | Amphar |
| 740313-3213 | 455 | 0.36 | | |
| 750108-3210 | 458 | 17.1 | 3 | |

Ex 31

C:\program\genesis\protocolab\307\NOV00V.004

LABSYSTEMS GENESIS 1.0.0.17
 Raw data filename: c:\program\genesis\protocolab\307\NOV00V.004
 Processed by Protocol: c:\program\genesis\protocolab\307\NOV00V.004
 Plotter layout file: c:\program\genesis\protocolab\307\NOV00V.004
 Reading type: Dual Wavelength
 Instrument version: Multiskan RC V.1.0
 Filter 1: 405nm
 Filter 2: 520nm
 Scan time: 00:00:10
 Interval between edge: 00:00:10
 Mix: YRS
 Mix RPM: High
 Mix ON period: 00:00:05
 Mix OFF period: 00:00:05
 Wavelengths:
 405 nm 520 nm

| Raw data values (calculated): | 1st | 2nd | 3rd | 4th | 5th | 6th | 7th | 8th | 9th | 10th | 11th | 12th |
|-------------------------------|-------|-------|-------|-------|-------|-------|-------|-------|-------|-------|-------|-------|
| A | 0.049 | 0.057 | 0.143 | 0.142 | 1.355 | 0.637 | 0.937 | 0.253 | 0.197 | 0.088 | 0.178 | 0.040 |
| B | 0.046 | 0.053 | 0.115 | 0.095 | 1.463 | 0.689 | 0.443 | 0.145 | 0.105 | 0.083 | 0.144 | 0.117 |
| C | 0.050 | 0.051 | 0.178 | 0.156 | 1.137 | 0.449 | 0.304 | 0.988 | 0.114 | 0.020 | 0.158 | 0.103 |
| D | 0.054 | 0.052 | 0.182 | 0.105 | 1.182 | 0.776 | 0.473 | 0.418 | 0.175 | 0.085 | 0.145 | 0.102 |
| E | 0.054 | 0.113 | 0.064 | 0.059 | 0.851 | 0.411 | 2.404 | 1.092 | 0.114 | 0.028 | 0.148 | 0.094 |
| F | 0.054 | 0.042 | 0.056 | 0.234 | 0.837 | 0.532 | 0.477 | 1.118 | 0.111 | 0.043 | 0.179 | 0.084 |
| G | 0.063 | 0.050 | 0.240 | 0.106 | 1.435 | 0.823 | 0.407 | 0.258 | 0.182 | 0.080 | 0.153 | 0.120 |
| H | 0.053 | 0.119 | 0.183 | 0.250 | 1.140 | 1.354 | 1.510 | 0.046 | 0.107 | 0.048 | 0.292 | 0.087 |

23 33 273 423 453 24 454
 11.4 285 26.5 2145 8.3

Amph 31

Comments:

Plate 3: 45 min

| Raw data values (calculated): | 1st | 2nd | 3rd | 4th | 5th | 6th | 7th | 8th | 9th | 10th | 11th | 12th |
|-------------------------------|-------|-------|-------|-------|-------|-------|-------|-------|-------|-------|-------|-------|
| A | 0.049 | 0.057 | 0.143 | 0.142 | 1.355 | 0.637 | 0.937 | 0.253 | 0.197 | 0.088 | 0.178 | 0.040 |
| B | 0.046 | 0.053 | 0.115 | 0.095 | 1.463 | 0.689 | 0.443 | 0.145 | 0.105 | 0.083 | 0.144 | 0.117 |
| C | 0.050 | 0.051 | 0.178 | 0.156 | 1.137 | 0.449 | 0.304 | 0.988 | 0.114 | 0.020 | 0.158 | 0.103 |
| D | 0.054 | 0.052 | 0.182 | 0.105 | 1.182 | 0.776 | 0.473 | 0.418 | 0.175 | 0.085 | 0.145 | 0.102 |
| E | 0.054 | 0.113 | 0.064 | 0.059 | 0.851 | 0.411 | 2.404 | 1.092 | 0.114 | 0.028 | 0.148 | 0.094 |
| F | 0.054 | 0.042 | 0.056 | 0.234 | 0.837 | 0.532 | 0.477 | 1.118 | 0.111 | 0.043 | 0.179 | 0.084 |
| G | 0.063 | 0.050 | 0.240 | 0.106 | 1.435 | 0.823 | 0.407 | 0.258 | 0.182 | 0.080 | 0.153 | 0.120 |
| H | 0.053 | 0.119 | 0.183 | 0.250 | 1.140 | 1.354 | 1.510 | 0.046 | 0.107 | 0.048 | 0.292 | 0.087 |

Coating: A: 10mg/ml
 B: 5mg/ml
 C: 2.5mg/ml
 D: 0.5mg/ml



cont. from page no:

cont. on page no:

Study

Project no.

Study no.

Test of Fd1 chain 1 clone 92 and
Fd1 chain 2, clone 29 with 6 cat sensitized
Allotopes from Allot study

6 farnes from the Allot study, no 29, 59, 277, 422,
434 and 454,
was diluted 3 times and 10 times resp. in
PBS pH 7.4. A μ -titer plate was coated with:

Plate 1

A 10 μ g/ml chain 1 \rightarrow (All horizontal wells)
B 5 μ g/ml "
C 2.5 μ g/ml "
D 1.25 μ g/ml "
E 10 μ g/ml chain 2 \rightarrow
F 5 "
G 2.5 "
H 1.25 "

Plate 2

A 5+5 μ g/ml chain 1+2 (a 1:1 mix of resp chain
B 2.5+2.5 "
C 1.25+1.25 " Chain 1+2 with 5 μ g/ml each)

On the vertical rows the patients were added

| | 1 | 2 | 3 | 4 | 5 | 6 | 7 | 8 | 9 | 10 | 11 | 12 |
|---------|--------------|--------------|--------------|--------------|--------------|--------------|--------------|--------------|--------------|--------------|--------------|--------------|
| | \downarrow | \downarrow | \downarrow | \downarrow | \downarrow | \downarrow | \downarrow | \downarrow | \downarrow | \downarrow | \downarrow | \downarrow |
| Patient | 29 | 59 | 277 | 422 | 434 | 454 | | | | | | |
| | 3x | 10x | 3x | 1x | | | | | | | | |

Read and understood by

Date

Signature

Date



cont. from page no.

cont. on page no.

Study

Project no.

Study no.

(cont.)

ELISA conditions

Coating of Fd 1:1 and Fd 1:2 ^(100 µl/well) over the weekend in +4°C. Serum 4 times was with "Tris" solution. Patients serum was added, 100 µl/well and incubated at +4°C o.n. Wash 4 times (Wallac ELISA-washer) with "Tris" solution and 100 µl Rabbit/human anti IgE dil 1000 times in "Vär buffer". Incubation 2 hours in RT on shaker, Wash 4x "Tris" solution and add 100 µl/well of Goat anti-rabbit-ALP conjugated (DAKO) for 1h. Wash 4 times and add substrate 3 tablets/15ml of Vär buffer.

The result was read in ELISA reader after 45 min at 405 ~~nm~~ nm

Result: 2.5 µg/ml seems to be an adequate coating concentration for both chain 1 and chain 2. Mixing of the two chains can be done with coating concentrations 2.5 + 2.5 µg/ml


Read and understood by

Signature

Date


Signature

Date

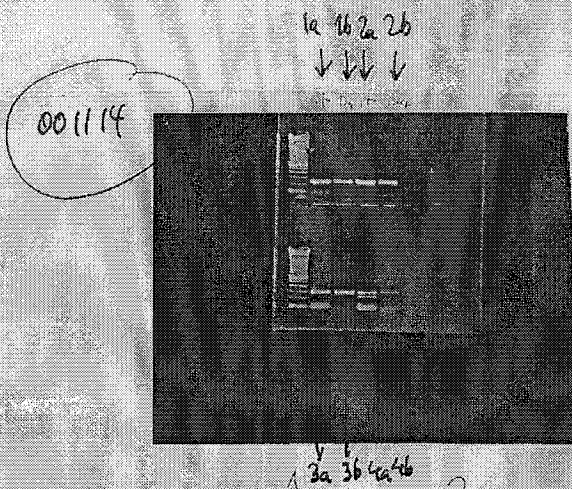

Putting chain 1 and 2 (1+2)
together (001121)

| | | |
|---------------------|--------|------------------|
| template 1a | 6µl | |
| template clone 29 | 1µl | |
| 1:10, 1:100, 1:1000 | 0.5µl | primer 176 |
| | 2µl | — — — 183 |
| | 2µl | dNTP |
| | 3µl | 10x buff. |
| | 1µl | Pfu |
| | 14.5µl | H ₂ O |
| | 30µl | |

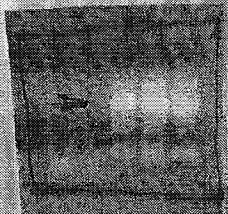

↑ ↑ ↑
1:10 1:100 1:1000

Result: One band
of ~500bp which
could be chain 1+2

↓
The bands are cut
out and purified
on Qiaquick

↓
ligated with "perfectly home
cloning kit". 10 colonies
are picked for miniprep
and possibly sequencing

①a 1µl template clone 11.42 (1:1000)
2µl primer 176
2µl — — — 174
2µl 10x buff
1µl Pfu
2µl dNTP (10mM)
10µl H₂O
20µl

①b Same as 1a, but
primer Tag polymerase

②a 1µl templ. (1:1000) clone 12
1µl templ. — — — clone 29
2µl primers. 176, 174, 183
overnight run ①a

②b Same ②a with Tag

③a 1µl template clone 2 (1:1000)
2µl 181
2µl 175
overnight run 1a

③b Same ③a with Tag

④a 1µl templ. 1
1µl templ. 2
2µl primer 186
overnight run ③a

④b Same ④a with Tag

AmpliTaq Gold
250 Units, 5U/µl
Store at -20 °C

Lot No. 00000000000000000000
A03912



cont. from page no:

cont. on page no:

Study

Linking of (chain 1 and chain 2 (seamless))
with PCR Fd1:1 and Fd1:2

Project no.

Study no.

Aim. The aim of this experiment is to join the two chains of Fd1 into one construct by PCR

Study outline The two sequenced chains of the major chains of cat Fd1 (chain 1, clone 42) and chain 2, clone 29) is joined with PCR in two steps as outlined below. In



Result (see opposite side)


Good bands of expected size was seen for both chain 1 and chain 2. But (2a) and (2b) as well as (4a) and (4b) did not work. I will continue by adding (1a) to ~~template~~ 29 and to PCR. (chain 1+2)



PCR: 5' (1+2)
1st templ. ch2 (1:1000)
2nd 186
3rd 184
→ did not work

PCR primers 180, 181

PCR: 5' (2+1)
1st templ. chain 1 (1:1000)
Chain 2 (1:1000)
1st 181
2nd 181
3rd 180
4th Pfu
5th 10x buff
6th 180
→ 5.30µl

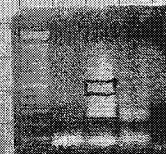

cut band

Read and understood by

Signature

Date

Signature



Date

PRISM











ABI106
Version 3.2

File 20
Lane 22

Seq 000
Points 19

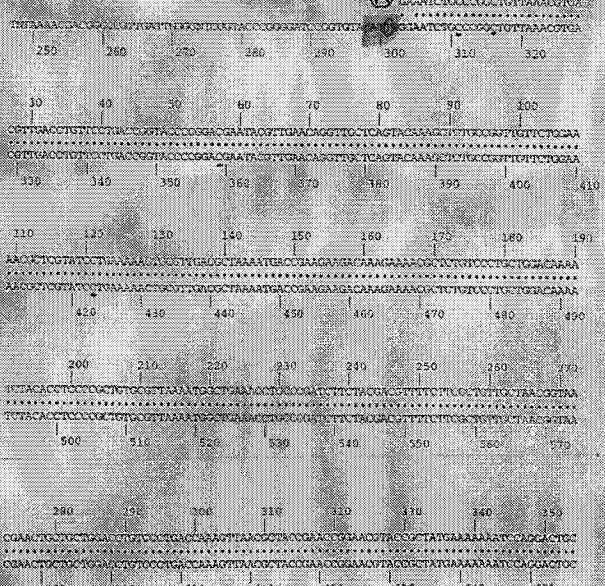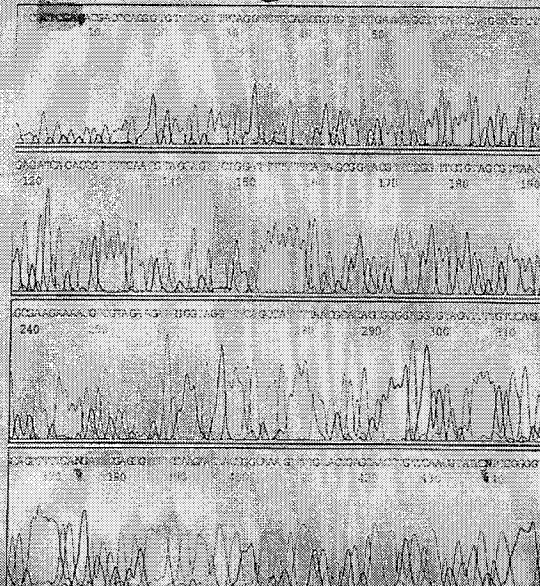

Ed1 (2+1) PT7 d1R 001205 (t+term.)

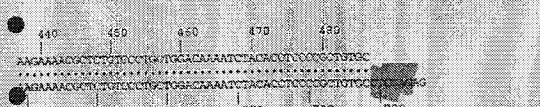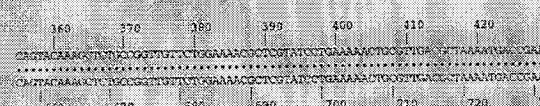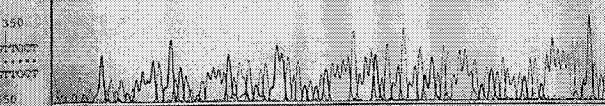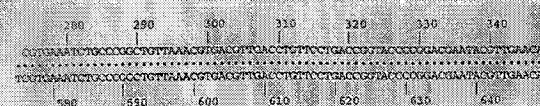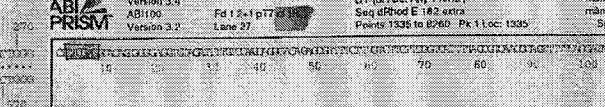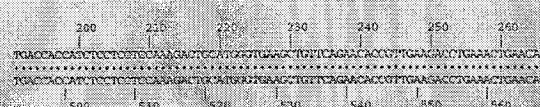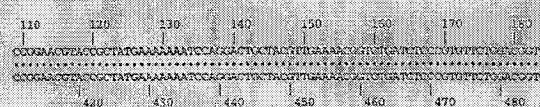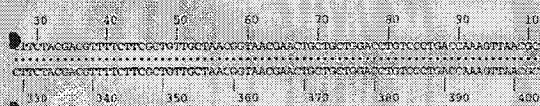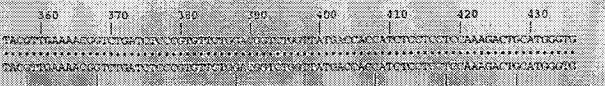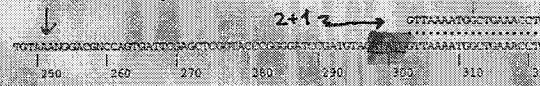


ABI
PRISM

Model 377
Version 3.4
ABI106
Version 3.2

File 20
Lane 22

Signal G112 A183 C111 T126
DT (dB) Set Any Primer
Seq diff 0.18 extra
Points 1335 to 8260 Pk 11 loc 1335

man 4 dec 8
man 4 dec 2
Spacing

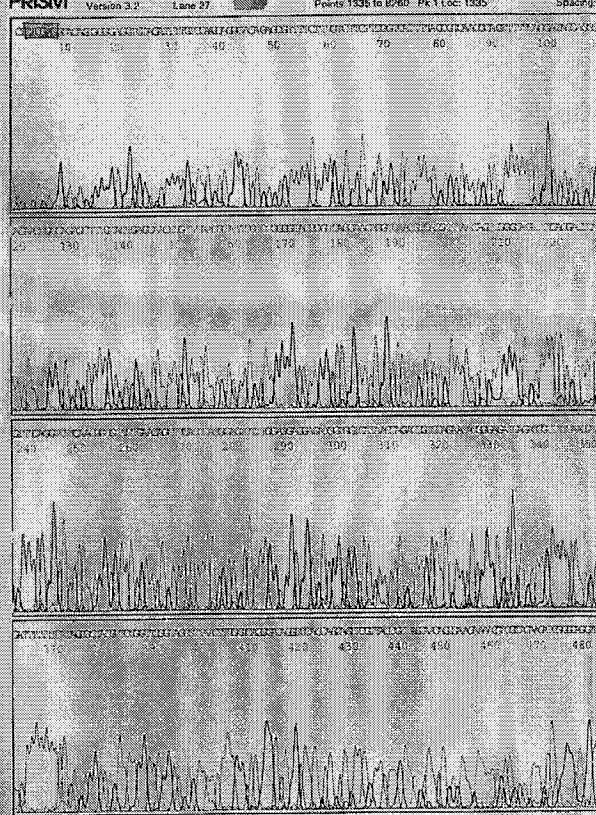



cont. from page no.

cont. on page no.

Study

Project no.

Study no.

Sequencing of 4 clones of each 1+2 and 2+1.

1+2 clones
↓
Clone 1 = ①
2 = ②
5 = ③ OK seq
9 = ④

2+1

Clone 1 = ⑤ OK seq red - prom
6 = ⑥ black - term.
7 = ⑦ OK seq (prom. by)
9 = ⑧

4.8 µl vector
1.2 µl primers
4 µl mix BD after dehydroamino
10 µl

+ 50 µl oil and

PCR p₀ Perkin Elmer

96°C 30"
50°C 15"
60°C 4" } 25 cycles

Unfortunately there was a scheduled power failure and the PCR-run was interrupted. Assume 7 more cycles which is done.

Samples are loaded on lanes 20 - 33 on ABI seq and named Fd 1 1+2 pT7 clone 4F

- " - 4R
- " - 5F
- " - 5R
etc.

Cloning of clone 5 (1+2) and clone 1 (2+1) with NotI and XhoI for ligation into pET 20b

20 µl plasmid miniprep (pT7 Blue)


2.4 µl 10x buff

0.3 µl BSA

0.7 µl NotI

0.7 µl XhoI

Incubate 37°C shake for 2h.

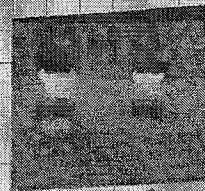

Diagonal point.

clone 1 (2+1)
clone 5 (1+2)

Read and understood by

Signature

Date

Signature

Date



cont. from page no.

cont. on page no.

Study

Ligation of Fd 1, clone 1 (2+1) and clone 5 (1+2)
into pET 20b+ and electroporation into BL21 plys

Project no.

Study no.


The fragments from Fd 1 (2+1) clone 1 and clone 5 (1+2)
were Qigen purified from 1% agarose gel (248931)
A cleaned (Nde/Xho) pET 20b+ vector was used
to ligate the fragments

Conclusions

12 μ l ~~vector~~ fragments
10 μ l vector
1.5 μ l 10mM ATP
1.8 μ l 10x ligase buffer
1.7 μ l T4 ligase
18 μ l

Ligate +16°C o.n.

The ligate mix was electroporated into 50 μ l BL21-
-plys electrocompetent cells. 1 μ l ligate mix
was added to thawed cells (on ice). Electroporation
according to standard protocol. Growth on SOC
medium for 60' 37°C shaker (300 rpm) and
plated on Amp/CAN plates. One colony ~~from~~
each plate was picked and grown on LB Amp
CAN medium, mini prepped (Qigen)
and 25 μ l of the (50 μ l) prep was
cut with Nde and Xho. Result ~ 800bp
Both clones contain the insert!!

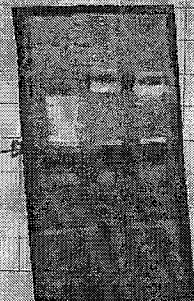

Read and understood by

Signature

Date

Signature

Date



cont. on page no:

Project no:

Study no:

Purification of Fd1 (1+2) clone 5 and
Fd1 (2+1) clone 1 over Ni^{2+} chelate thi-Trap

1 liter of Fd1 (1+2) and (2+1) ^{in PET 20b in BL-2 (Phys)} resp. LB-medium (CAM, Amp)
was grown to OD 0.6 (600nm) and induced with
0.4 mM IPTG. (see 248932)

Purification according to protocol. Both proteins
were expressed as inclusion bodies and
purified accordingly ~~and~~ Purification on FPLC
as follows. After adsorption onto thi-column in
6M Urea and wash also with 6M Urea the
column is stuck to FPLC

Program:

0 conc % B 0

0 ml/min 5.0 ml/min

0 0.25 ml/min

0 port set 6.1

20 conc % B 0

80 conc % B 100

100 conc % B 100

125 conc % B 0

125 port set 6.0

A = 6 M Urea

B = 20 mM Imidazole

C = 500 mM Imidazole



Read and understood by

Signature

Date

Signature

Date